



Explorer endoscopy face masks

For analgesic gas administration during endoscopy procedures



Anaesthesia ▪ Anaesthetic Face Masks

Explorer endoscopy face masks

For analgesic gas administration during endoscopy procedures

A range of full face PVC masks for analgesic gas administration during endoscopy procedures, which are single use to prevent the possibility of cross infection.

The Intersurgical Explorer endoscopy masks are designed to allow easy insertion of an endoscope probe into the nose or mouth whilst administering an analgesic gas. There are three one-way valves positioned on the mask (the size 1 mask has only one valve) which provide easy access from either side of the patient.

Principally used for bronchoscopy procedures, the analgesic gas is administered via the central 22mm port. These masks could also be used on patients undergoing a transoesophageal echocardiogram (T.O.E.) procedure.

▶ Video available at www.intersurgical.com



Colour-coded hook ring

For ease of identification, and integration with other Intersurgical airway products

Self-sealing valves

To facilitate easy access

Anatomical shape


For a precise fit

Soft cushion

For patient comfort and air tight fit

4 sizes available

From adult to infant

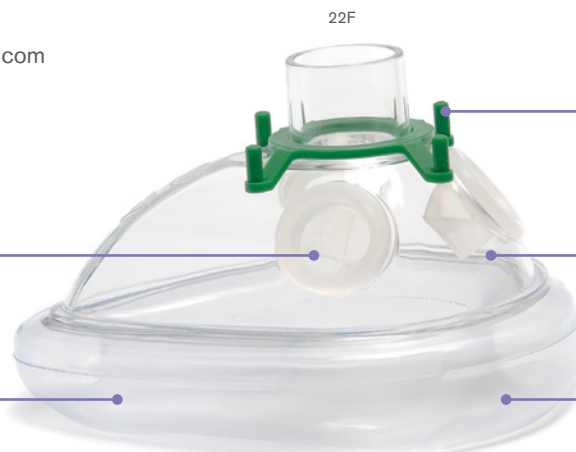

Code	Description	Box Qty.
1594000	Explorer, endoscopy face mask, size 4, adult, 22F	40
1593000	Explorer, endoscopy face mask, size 3, small adult, 22F	50
1592000	Explorer, endoscopy face mask, size 2, paediatric, 22F	35
1591000	Explorer, endoscopy face mask, size 1, infant, 15M	40

Make an enquiry



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

Please think before you print
 Save energy and paper.
 If you must print this information sheet please print it double sided.

IS5.18 • Issue 5 07.19