



Intersurgical EcoLite™ adult high concentration oxygen mask



Oxygen and Aerosol Therapy ▪ Oxygen Therapy



Intersurgical EcoLite™ adult high concentration oxygen mask

Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials in the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patients face.

Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask body and seal, reducing the environmental impact.

Intersurgical EcoLite

[Find out more](#)

Eco range

[Find out more](#)



Incurved nose seal

Conforms to different nose shapes designed to minimise the risk of oxygen entering patient's eyes

No metal nose clip

MRI compatible

Clear mask

permits visualisation of the patient's nose and lips

Two integrated safety valves

Allows patient to breathe room air in the event oxygen supply failure

Soft face seals



Increased patient comfort

Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears

A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes

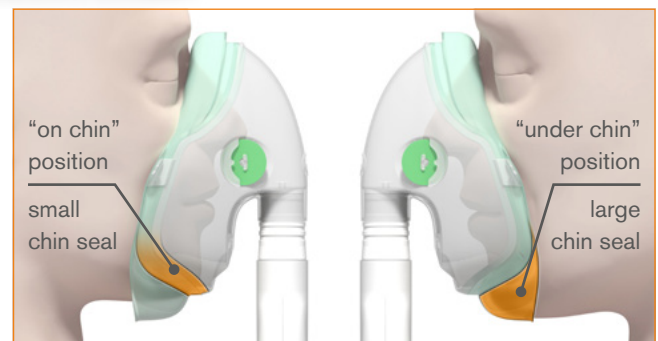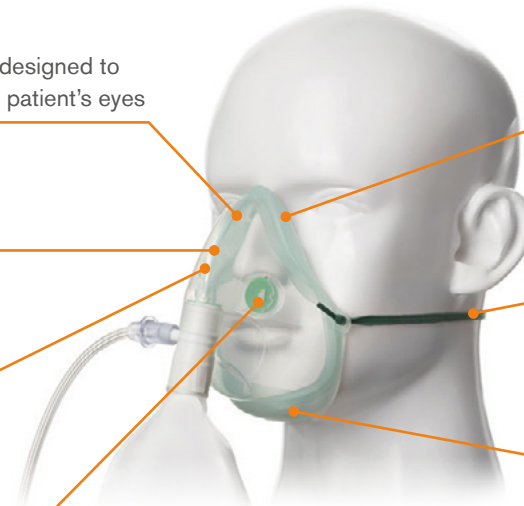

Code	Description	Tube length	Box Qty.
1181015	Intersurgical EcoLite, adult, high concentration oxygen mask with tube	2.1m	24

[Make an enquiry](#) [Watch the video](#)

Lower environmental impact product. Note: reservoir bag and tube contain PVC.

IS10.17 • Issue 11 03.24



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

Please think before you print
Save energy and paper.
If you must print this information sheet please print it double sided.